


la belle BC 211 P

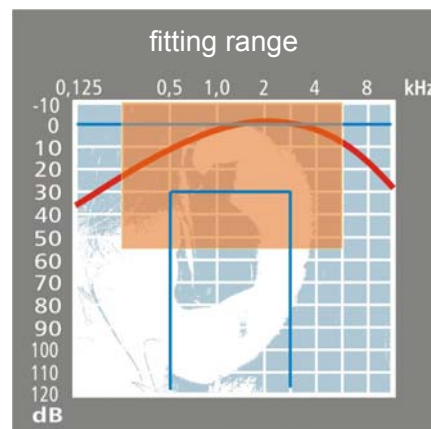


la belle BC 211 P the bone conductor glasses

- BTE-module for eyeglasses
- module cold bendable
- digital hearing-system
- 3 programmes selected by push button
- 2-channel WDRC
- Dynamic Contrast Detection
- 12-band gain adjustment
- volume control VC
- telecoil
- MPO
- adaptive feedback cancellation (until -17dB)
- external microphone on spectacle arm
- receiver optionally above or below (normally below)
- programme tones and low battery tones
- battery type 312
- colours: anthracite-black
- colour of the chamber cap: silver, gold, anthracite
- feeler pin in support of finding the appropriate glasses

		DIN IEC 60118-9:1987
FOG for OSPL90 (1600 Hz)	(db)	118
NAW (1600 Hz)	(dB)	53
Frequency range	(Hz)	250-5500
Operating current	(mA)	1,2
Equivalent input noise	(dB)	22

please turn over

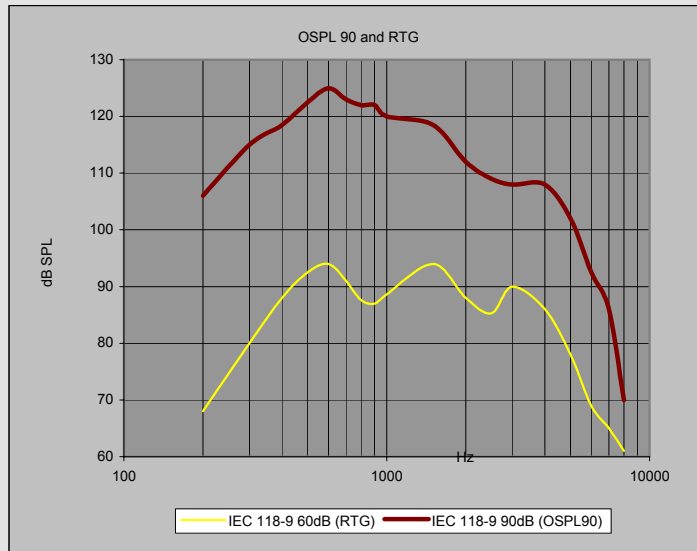

Programming:

with Bruckhoff fitting software and Siemens programming cable (CONNEXX cable grey)

Technical change without prior notice